



Plasma clean
Improved wettability and adhesion
Functional nano-coatings

Plasma Surface Treatment

HPT-300 Benchtop plasma treater



HPT-300 Benchtop Plasma Treater

The HPT-300 is a microprocessor controlled benchtop plasma treatment system which is ideally suited to surface activation, cleaning and modification of a wide range of materials including polymers, metals, glass and ceramics.


Available in single or dual gas inlet versions and with on-board gas mixing manifold, the HPT-300 is able to handle a wide range of gases for optimized treatments, including air, oxygen, hydrogen, argon, nitrogen and many others.

An optional vapor delivery inlet extends the use to liquid precursors and a corrosion resistant version expands the choice even further to address specific material treatments including;

- Plasma cleaning
- Plasma surface activation to improve adhesion
- Functional plasma coatings
- Plasma etching
- PDMS & microfluidic devices
- PEEK & other engineering polymers
- PTFE
- Metals
- Ceramics
- Glass & optical devices



Markets & Applications

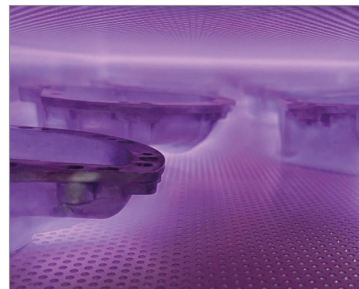

**Aerospace &
Automotive**

Composites



Plasma Environment



Plasma Environment

The HPT-300 features a 190mmW x 310mmL plasma process chamber in stainless steel with vacuum compatible materials throughout. Our proprietary, high stability HPS plasma generator is continuously variable over the entire 0-250W power range rather than being limited to discreet levels, delivering much finer control when processing delicate materials.




Process Control

The 5.7 inch color touchscreen provides a rich, user-friendly interface. Variables such as gas flow rate, pressure, power level and plasma processing time can be freely set and then stored to produce a fully interlocked process cycle from a single keypress. A handy status display and end of process audible alarm informs the user of every step in the process.




Repeatable & Reliable

With precision digital mass flow controllers and integrated pressure gauge, the HPT-300 delivers unmatched reliability and repeatability by removing common errors in gas flow and gas type settings which will be familiar to users of equivalent equipment that utilize manual needle valves. Convenient recipe store software, a unique feature, allows fixed repeatability.

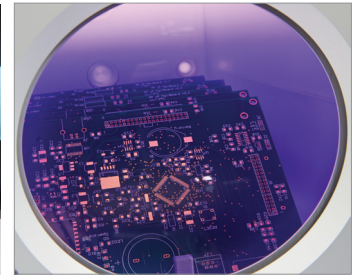


Versatile


Recipe store allows users to select and store up to four recipes. Each recipe contains unique settings for power level, plasma process time and pressure. The base model HPT-300 has a single gas inlet and optional second gas or vapor delivery inlet. The unit is prepared so that either option can be added at a later time if required, ensuring that future requirements can be accommodated without expensive reconfiguration.

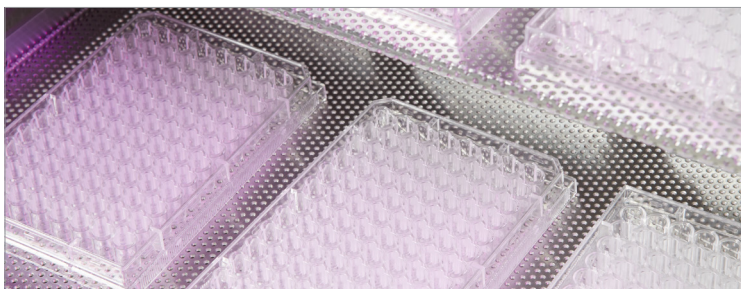

Medical Plastics

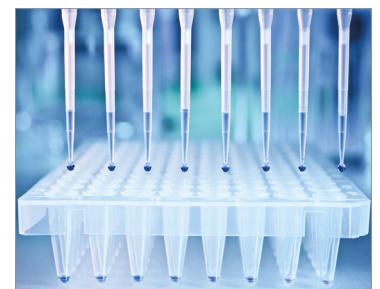

Microfluidics

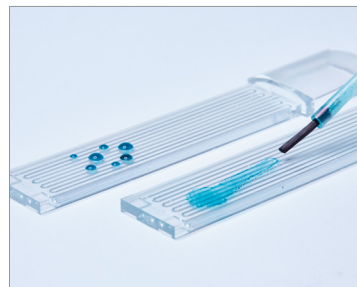

Some of Our Clients



"Henniker provided visible results from the outset and confirming that we made the right decision in choosing a local UK manufacturer."

Queen's University Belfast



"We obtained quality results with their unit within minutes of setup & consistent results thereafter. The support they have provided has been rapid and helpful."

Making Lab, Francis Crick Institute



"We are very impressed with the ease of use and reliability of our plasma unit and were producing results within minutes of setting it up."

Warwick University



"Our collaborative work with the team at Henniker was a very positive experience and one that we look forward to developing further."

TWI



TAGHeuer



KONGSBERG
AUTOMOTIVE



i n v e n t




**UNIVERSITY OF
CAMBRIDGE**



**Plasma
Cleaning**



**In-Line
Integration**

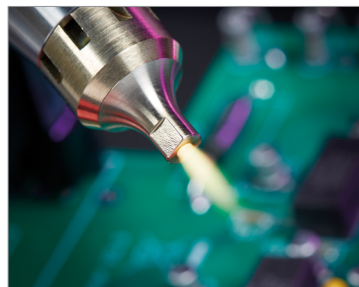

PCBs



HPT-300 Specifications


BASE MODEL		OPTIONS
ENCLOSURE		
Dimensions	W 533mm x H 466mm x L 615mm (+50mm on rear for cables)	
Weight	~85lbs	
CHAMBER		
Material	Stainless Steel	
Form	Rectangular	
Dimensions	200mm W x 200mm H 320mm L	
REMOVABLE PARTS CARRIER		
Material	Aluminum	stainless steel
Form	Flat tray	multi-level shelf/electrode
Dimensions	190mm W x 310mm L	others to suit application
PLASMA POWER SUPPLY		
Power	0-250W, continuously variable output	
Frequency	40 kHz	
PROCESS CONTROL		
Interface	5.7" Color TFT with recipe store	
Gas channels	x1 MFC	x2 MFC or x1 MFC and x1 vapor inlet
Vent inlet	x1	soft ventilation option
Connections	1/4" compression	
Process timer	1sec – 99.59min	
Recipe Store	Stores up to 4 x recipes with individual parameters	
Pressure gauge	Pirani sensor	
Vacuum pump	6 to 15 m3/hr pumping speed	
Vacuum pump options	2-stage rotary pump (air/inert gas), PFPE rotary pump (oxygen compatible), dry pumps. All pumps include exhaust filter and connections	
SERVICES		
Electrical	90-120 VAC, 50-60Hz, 1200 VA (including pump), fused 10 A T	
Power cord	6ft	
Compliance	CE – UKCA - ROHS - WEEE	

Benefits


- compact benchtop unit
- user friendly TFT interface
- recipe store
- fast treatment time
- precise & repeatable
- no hazardous emissions

Typical Process Results

- Material : PEEK
- Oxygen plasma
- 60 seconds

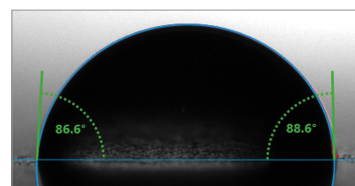

Contact Angle Before Treatment

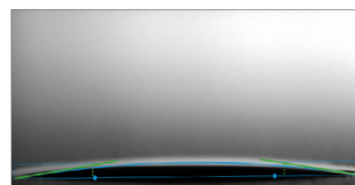


Contact Angle After Treatment


Henniker strive for continuous improvement and specifications are subject to change without notice



Products & Services

- benchtop systems
- high throughput systems
- atmospheric plasma
- robot systems
- surface test & analysis
- process development

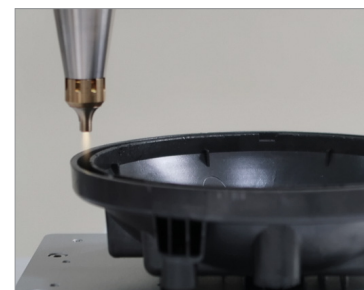

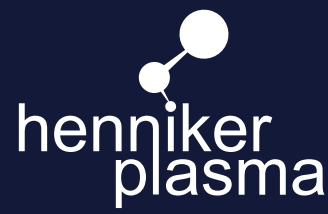

About Henniker

Henniker Plasma are an international leader in the design, development and manufacture of plasma surface treatment systems & advanced plasma processes.

Our products are installed worldwide and trusted to deliver consistent, reliable results in both leading research institutes and in critical manufacturing steps.

We are experts in plasma technology and surface science. We are trusted partners, valued for our courtesy, professionalism and dedication to delivering the correct solution for our clients.

Services

Contract plasma treatment

Our technical staff will be happy to discuss contract treatments, from small one-off batches to regular, large throughput requirements.

Proof of concept treatment

Let's discuss your application and then we will provide a quick, no-nonsense feasibility study.

Surface testing laboratory

With a comprehensive suite of surface analysis equipment, we are able to conduct a wide range of surface property tests, both before and after plasma treatment, in order to provide you with the whole picture.

After sales support

We are proud of our reputation for being approachable, thorough and easy to work with.

"Henniker's after sales support is first class. They have always been extremely responsive if we have ever had need to call on them."

Steve Rackham, Teledyne

Rental plasma systems

We carry a wide range of our standard equipment in stock and available for short or long term hire. This is particularly useful for in-house proof of concept trials or to satisfy short term contract work.

"The low risk option of hiring a plasma unit for evaluation was a key reason that we chose to work with Henniker and one that enabled us to proceed with confidence."

Dr. Chris Nicklin, Reinnervate

Method development

We have invested significantly in laboratory facilities to assess, test and investigate all aspects of plasma surface modification on a wide range of materials. Coupled with extensive in-house and real-world knowledge, we can usually deliver a tailored treatment quickly and efficiently to suit your individual product or production needs.

"The technical team at Henniker are very knowledgeable and supportive and always approachable. I have found it a pleasure to work with them."

Simon Baxter, BAE Systems, AI

Henniker Plasma

3 Berkeley Court
Manor Park
Runcorn WA7 1TQ
United Kingdom
Tel: +44 (0)1925 830 771
Web: www.plasmatreatment.co.uk



Represented by:

**PRINCETON
SCIENTIFIC**
CORPORATION

P.O. Box 148 • Easton, PA 18044
Tel: (609) 924-3011 • Fax: (609) 924-3018
www.PrincetonScientific.com
Email: info@princetonscientific.com